

SINGLE USE GLOVE, 100% NITRILE, POWDER FREE, 240MM LONG, 0.10MM THICK

High performance, powder-free, nitrile single use examination glove. Significantly stronger alternative to latex, dual certified (Industrial and Medical) and food safe.

FEATURES

- › High quality nitrile, soft and durable
- › Palm thickness 0.10mm (4 mil), weight: 5g size M
- › High resistance to penetration and projections of many chemicals
- › Textured finish on fingertips for improved grip and touch
- › Powder-free, latex free, non-sterile, AQL 1.5
- › Chlorinated inside to offer more comfort and reduce friction when donning
- › Antistatic properties due to inherent antistatic properties of nitrile
- › Approved for food handling
- › Dual labelling: PPE and medical
- › Medical device compliant (MDD 93/42/EEC), conforms to EN455 (all parts)
- › Tested against chemotherapy drugs according to ASTM D6978-05
- › Compliant with ASTM D6319-2019 and ASTM D5151-2019



SUITABLE FOR


Typical Industries

- Aerospace
- Chemical
- Electronics
- Food
- Pharmaceutical
- Manufacturing
- Oil and Gas

Suitable Applications

- Chemical Handling
- Food Handling
- Handling Oily Components
- Laboratory Analysis
- Light Assembly
- Paint Spraying
- Examination

CERTIFICATION

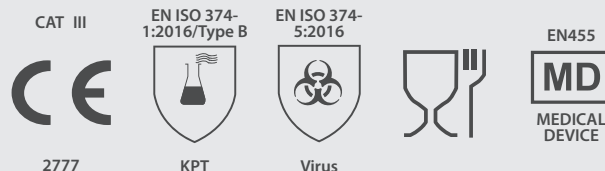

See overleaf for explanation

PRODUCT INFORMATION

MATERIALS	LINER:	Nitrile
	COATING:	N/A
COLOUR	Blue	
LENGTH (mm)	240	
CUFF STYLE	Beaded	

ORDERING DETAILS

SIZE	CODE	PACKAGING
6/XS	SKG00424DB	100 singles per box (90 singles on size XXL)
7/S	SKG00424DD	
8/M	SKG00424DF	
9/L	SKG00424DH	
10/XL	SKG00424DJ	
11/XXL	SKG00424DL	10 boxes per case

RECOMMENDATIONS FOR USE

- **USE:** Single use glove. Not suitable for mechanical, thermal, heat, electrical, heavy chemical protection. Do not use near moving machines if there is a risk of entanglement
- **STORAGE:** Store in dry conditions in the original packaging and away from direct sunlight
- **CLEANING:** To clean, wipe with a damp cloth. Note: The performance characteristics of worn and laundered gloves may differ from the results shown. Inspect the gloves to ensure no damage is present.
- **LIFETIME:** Service life depends on the glove application and therefore cannot be specified. It is the responsibility of user to ensure the glove is suitable for its intended use

CERTIFICATION LEGENDS



RESISTANCE TO CHEMICAL PERMEATION - EN ISO 374:2016

CODE	CHEMICAL	CODE	CHEMICAL	TYPE OF GLOVES	BREAKTHROUGH TIME
A	Methanol	J	n-Heptane	A	≥30 min for at least 6 chemicals
B	Acetone	K	Sodium hydroxide 40%		
C	Acetonitrile	L	Sulphuric acid 96%		
D	Dichloromethane	M	65% Nitric acid	B	≥30 min for at least 3 chemicals
E	Carbon Disulfide	N	99% Acetic acid		
F	Toluene	O	25% Ammonium hydroxide		
G	Diethylamine	P	30% Hydrogen peroxide	C	≥10 min for at least 1 chemical
H	Tetrahydrofuran	S	40% Hydrofluoric acid		
I	Ethyl acetate	T	37% Formaldehyde		

EN455-1:	Freedom from holes
EN455-2:	Physical properties
EN455-3:	Biological evaluation
EN455-4:	Shelf Life



PROTECTION AGAINST MICRO-ORGANISMS EN 374-5

VIRUS = Glove has passed ISO 16604:2004 (method B)



FOOD CONTACT DIRECTIVE EC 10/2011

EUROPE

www.globusgroup.com

E: sales@globus.co.uk
T: +44 (0)161 877 4747
F: +44 (0)161 877 4746

MIDDLE EAST AND AFRICA

www.globusgroup.com/gcc

E: gcc@globusgroup.com
T: +971 4 882 9962
F: +971 4 882 9963

